

Arterion™ offers multiple configurations for maximum configuration flexibility



Integrated Pedestal
The integrated pedestal has been completely re-designed for function and style



Table Mount
The small power supply allows mounting flexibility – under the table, in control room, in store room



Freestanding Pedestal
Ideal for moving the injector head on and off table when not needed



Overhead Counterpoise System
The ceiling mount option allows the injector head to be moved easily for injections

Option for all systems:

Heat Maintainer

Keep contrast temperature constant with Arterion's clear heat maintainer



Accessorize Arterion to meet your unique clinical needs

Dual Display Add a second display for maximum operational flexibility between patient room and control room

SYSTEM TECHNICAL SPECIFICATIONS:

Flow Rate:	0.1-45.0 ml/s in 0.1 ml increments 0.1-59.9 ml/m in 0.1 ml increments
Volume:	1-150 ml in 1 ml increments
Pressure Limit (150 ml syringe):	100-1200 psi in 1 psi increments 689-8273 kPa in 1 kPa increments
Rise Time:	0.0-9.9 seconds in 0.1 increments
Delay Time:	0.0-99.9 seconds in 0.1 increments
Fill Speed:	1-20 ml/s
Fill Volume:	1-150 ml
Syringe Size:	150 ml
Protocol Memory:	40 Protocols
Injection Memory History:	50 Injections

Contact your local MEDRAD Representative for more information on a Mark 7 Arterion to meet your unique needs.



The Mark 7 Arterion Injection System is marketed for the administration of intravascular contrast media and common flushing agents during angiographic procedures. See Product Information for Use for specific and complete prescribing information. Unless otherwise indicated, all trademarks are owned by MEDRAD, INC. or licensed for its use. MEDRAD Interventional is a federally registered trademark of MEDRAD, INC. U.S.A. ©2011 MEDRAD, INC. All Rights Reserved.

MEDRAD®
Performance. For life.™

MEDRAD, INC.
One Medrad Drive
Indianola, PA 15051-0780 USA
Ph: +1 (412) 767-2400
Other: +1 (800) 633-7231
Fax: +1 (412) 767-4120
www.medrad.com

MEDRAD Europe B.V.
P.O. Box 205
6190 AE Beek
The Netherlands
Ph: +31 (0) 43-3585600

A Business of Bayer HealthCare

2094-001 Rev 3 11/11



MARK 7
ARTERION™
Injection System

We deliver contrast.
You deliver exceptional patient care.

MEDRAD
Interventional®



ergonomic handle
for easier maneuverability

unique front
load syringe
for simple insertion
and clean removal

16% lighter
injector head*

32% lighter
pedestal unit*
Light and easy to move

intuitive
interface
for simple
injector control

smooth
arc design
for extended reach

20% smaller
footprint*

Mark 7 Arterion has stylish ergonomics that are designed for function; the pedestal system with a smaller footprint is lighter than MEDRAD's previous version, which increases mobility around your busy lab. And the lighter head now has a handle to make it easier to position for injection.

The Mark 7 Arterion™ Injection System is the latest in MEDRAD Interventional's "Mark" series of angiographic injectors.


The Mark 7 Arterion takes advantage of the latest technologies, making it light, maneuverable and easy to use. Less time positioning and setting up the Arterion means more time with the patient. The clearly visible and intuitive user interface guides you through proper set-up, and highlights the information you need to perform injections confidently.

Unique to the market, the front load system simplifies set-up and makes for a cleaner tear down. The Mark 7 Arterion syringe provides a clear view of the contrast. With a smaller base, lighter head and a flexible articulating arm, the injector can go where it needs to go – even around corners.

And of course you can continue to count on what MEDRAD has always delivered: Reliable equipment with an extensive team of highly trained people ready to provide exceptional service.



Made from a clear material, the Mark 7 Arterion syringe (Catalog ART 700 SYR) allows you to easily view the inside of the syringe for monitoring the purging of air. And MEDRAD's famous fluid dots are still there to help—round for fluid, oval for air.

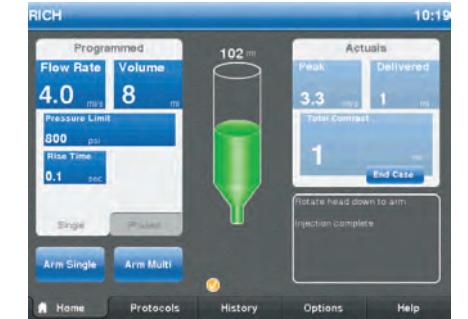

DISARMED STATE
The bright and colorful user interface is designed to highlight the information you need.



ARMED STATE
The highlighted armed state also lets you know when the system is ready to inject. Our modified history and protocol screens will allow you to easily access the amount of contrast delivered to the patient and store and recall protocols.



*as compared to ProVis®