

Administered 95.9 ml of Omnipaque 350 mg/ml via 18 g with P3T Cardiac protocol



Injection Start: 24 May 2010 11:31 am
 Accession Number: 555123888
 Patient ID: 000333222



Image Created: 24 May 2010 11:48 am
 Patient Name: Mary Jones
 Date of Birth: 01 March 1950

P3T Cardiac

Programmed			Actual						
	ml/s	ml		ml/s	ml				
1	B	5.7	20	1	B	5.2	19.4		
2	Hold			2	Hold				
3	A	5.7	20	3	A	5.3	19.4		
4	B	5.7	40	4	B	5.3	39.9		
5	A	5.7	68	5	A	5.5	67.7		
6	A	20%	5.7	46	6	A	20%	5.5	44.3
7	B	5.7	30	7	B	5.6	29.9		

Total Contrast (A): 97 ml
 Total Saline (B): 127 ml
 Delay: Scan Delay: 20 s
 Pressure Limit: 325 psi

Total Contrast (A): 95.9 ml
 Total Saline (B): 124.7 ml

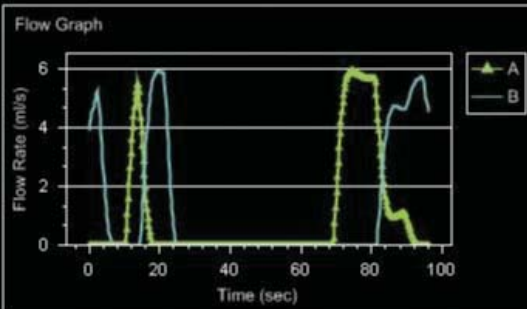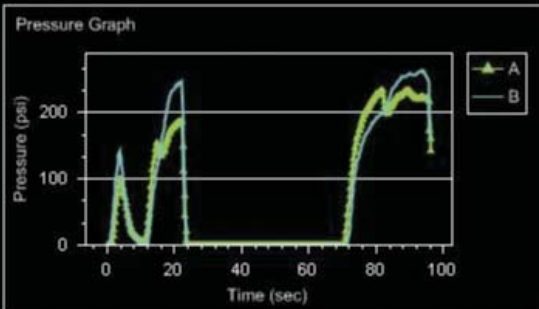

Injection	Contrast	Saline
Peak Pressure:	238 psi	265 psi
Peak Flow Rate:	6.0 ml/s	6.0 ml/s
Injection Completed:	YES	
Transient Events:	NO	

Fluid Usage	Contrast	Saline
Loaded:	106.3 ml	162.3 ml
Used:	95.9 ml	124.7 ml
Remaining:	10.4 ml	37.6 ml
Total Fluid:	220.6 ml	
Total Iodine:	33.6 gI	

Contrast
 Brand: Omnipaque
 Concentration: 350
 Lot Number: 10827593
 Expiration Date: 31 December 2012
 Vial Volume: 500

Procedure
 Tech ID: slw
 Injection Site: RAC
 Catheter Gauge: 18
 Injector Name: Injector

Notes
 --



Not actual patient information. For demonstration purposes only.

Administered 69.9 ml of Omnipaque 350 mg/ml via 22 g with P3T PA protocol



Injection Start: 24 May 2010 10:37 am
 Accession Number: 668772321
 Patient ID: 000333222

Image Created: 24 May 2010 10:57 am
 Patient Name: Mary Jones
 Date of Birth: 01 March 1950



P3T PA

Programmed		ml/s	ml
1	B	4.8	20
2	Hold		
3	A	4.8	57
4	A 20%	4.8	67
5	B	4.8	30

Actual		ml/s	ml
1	B	4.4	19.4
2	Hold		
3	A	4.5	56.9
4	A 20%	3.5	65.3
5	B	4.3	30.4

Total Contrast (A): 70 ml
 Total Saline (B): 104 ml
 Delay: None
 Pressure Limit: 325 psi

Total Contrast (A): 69.9 ml
 Total Saline (B): 102.0 ml

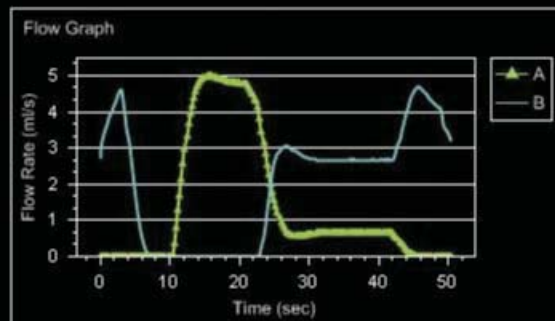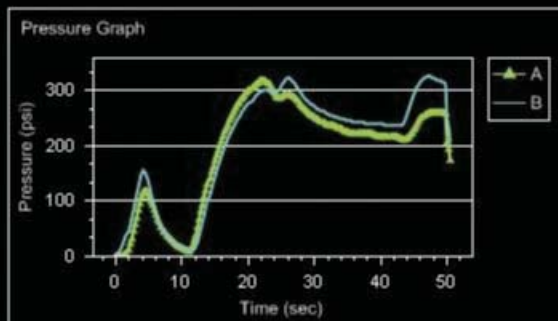

Injection	Contrast	Saline
-----------	----------	--------

Peak Pressure: 326 psi 333 psi
 Peak Flow Rate: 5.1 ml/s 5.0 ml/s

Fluid Usage	Contrast	Saline
-------------	----------	--------

Loaded: 72.3 ml 111.5 ml
 Used: 69.9 ml 102.0 ml
 Remaining: 2.4 ml 9.5 ml
 Total Fluid: 171.9 ml
 Total Iodine: 24.5 gI

Injection Completed: YES
Transient Events: Flow Rate Limiting

Contrast

Brand: Omnipaque
 Concentration: 350
 Lot Number: 10827593
 Expiration Date: 31 December 2012
 Vial Volume: 500

Procedure

Tech ID: ljack
 Injection Site: RAC
 Catheter Gauge: 22
 Injector Name: Injector

Notes

SOB



Not actual patient information. For demonstration purposes only.