

New Disposables Offering Continuum MR Infusion System

The new Continuum microbore disposables simplify MRI infusion workflow and make it easier for you to provide MRI access to your medication-dependent patients. The disposables have been designed to extend your capabilities to:

- Increase the types and volume of applications performed
- Provide more immediate patient care
- Lower procedure cost by reducing required pharmaceutical priming volume
- Enable backup medication source; visually detect air through the clear y-connectors.

Primary Microbore

MEDRAD part number 3015153 – Priming Volume: ~7.5 mL


- Current MIK200A set design with Microbore and Tru-swab
- Use case examples
 - Direct connect to critical care unit patients with ongoing IV therapy
 - Emergency Care patient



Secondary Microbore

MEDRAD part number 3015155 – Priming Volume: ~3.0 mL

- Current MIL200B set design with Microbore and Tru-swab
- Use case examples
 - Piggy back multi medication into primary line
 - Direct connect to existing meds

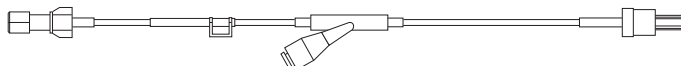

Primary with VSA plus Stopcock

MEDRAD part number 3015156 – Priming Volume: ~7.5 mL


- New Primary set with vented syringe adapter and stopcock
- Use case examples
 - Syringe source application with option for syringe backup or flush
 - Sedation with backup syringe source of same med
 - Medication delivery using one pump configuration with saline flush



Primary plus VSA

MEDRAD part number 3015157 – Priming Volume: ~7.5 mL

- New Primary set with vented syringe adapter
- Use case examples
 - Single syringe source application

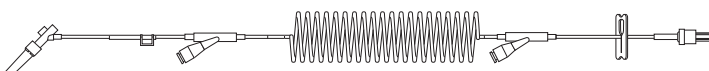

Primary plus Spike

MEDRAD part number 3015158 – Priming Volume: ~7.5 mL



- New Primary set with Vented Spike
- Use case examples
 - New bag source applications started with Continuum
 - Emergency/Trauma care opportunity
 - Start new medications in the MR suite or on the floor prior to transfer



SideStep

MEDRAD part number 3015159 – Priming Volume: ~4.5 mL


- New set attached to bypass the floor disposable key set without disconnecting the patient
- Use case examples
 - Any current peristaltic pump disposable that has a y-connector above and below the key / cassette
 - Typically ICU patients on infusion therapy

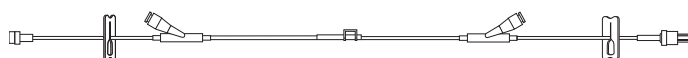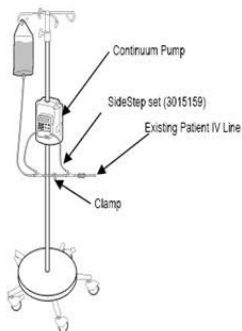

Secondary plus Spike

MEDRAD part number 3015160 – Priming Volume: ~3.0 mL


- New Secondary set with Vented Spike
- Use case examples
 - Capability to piggy back onto primary sets
 - Multi bag application common use Saline, D5W

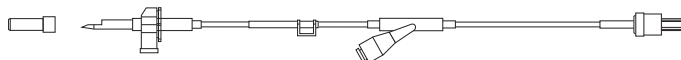

Secondary plus VSA

MEDRAD part number 3015161 – Priming Volume: ~3.0 mL

- New secondary set with Vented Syringe Adapter
- Use case examples
 - Piggy back to Primary set with stopcock/VSA
 - Sedation cases using syringe source applications
 - Piggy back syringe with medication to primary
 - Medication delivery using a two pump configuration

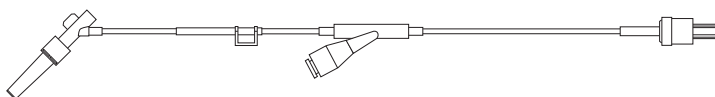

**For additional information on MEDRAD's Continuum Disposable Portfolio
or for ordering information, contact MEDRAD Customer Service at**

1-800-633-7231