



Continuum Wireless Infusion Pump Upgrade Program



Program Details:

For a limited time only, MEDRAD is offering a special promotion giving customers the opportunity to upgrade to our latest MRI infusion technology at a discounted price. The offer includes:

- Hardware and software feature enhancements
- A brand new pump with "wireless ready" charging bracket and a 12 month warranty; this "wireless ready" feature positions your site to purchase the optional Wireless Remote Display at any time
- An additional 5% discount on the Continuum Wireless Remote Display when it is purchased during the promotional period (this discount is in addition to your contracted price)

This special upgrade program ends August 31, 2010.

Beginning September 1, 2010, customers wishing to take advantage of these new hardware or software features, or those wishing to purchase a wireless pump and Remote Display will be required to purchase a new infusion system at their contracted system price.

Hardware and Software Enhancements offer Safety and Workflow Benefits:

Hardware Enhancements:

- New Bracket Design
 - Allows for In Bracket Charging, enabling you to perform stack charging of multiple pumps with one power supply
 - Eliminates the need to remove the pump from the bracket to charge the pump
 - Secure locking design reduces the potential for staff to remove the pump inside the scan room
- MR Compatible power supply allows you to charge and operate the pump while performing procedures inside the scan room


Software Enhancements

- KVO setting range expanded to 0-5ml/hr. Zero turns KVO off, which is a benefit of particular importance for sensitive cases involving smaller patients.
- Infusion Near End feature
- Dose/Rate Unit Lock: Safety feature that allows user to lock a single medicine for a given pump
- Bolus Review Screen: Safety enhancement that allows user to review bolus dose before delivery

In addition to the hardware and software features described, upgrading to the Continuum Wireless Ready Pump positions your site to take advantage of all of the safety and workflow benefits of the optional Continuum Wireless Remote Accessory.

The Remote Display offers:

- A consistent and reliable connection from the Infusion Pump inside the scan room to the Remote Display positioned in the control room
- Easy to follow screens which streamline throughput and provide user confidence in infusion management outside the scan room

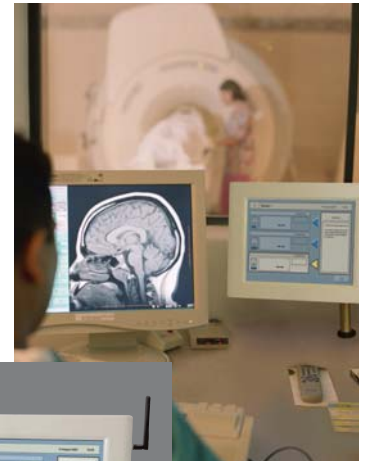

Continuum Wireless Remote Benefits:

- Eliminates inefficient trips into the scan room by allowing command and control at the Remote
- User can titrate flow rates, program and deliver bolus, and start or stop an infusion at the Remote
- Provides immediate visibility to audible and visual alarms, allowing the user to react more quickly to patient needs
- Flexibility to monitor up to nine pre-configured stands
- Ability to program drug name into the remote, providing improved visibility to all staff about what drugs are being administered

These enhanced benefits of Continuum Wireless are all in addition to the core features that the Continuum MRI Infusion System has built its name on:

- Three Infusion Modes: Continuous Flow, Weight Based Dosing, Non-Weight Based Dosing
- Dose units: mg, mcg, U, mU with time units of min. or hr.
- Bolus mode with safe and accurate delivery interface
- Flow rates: 0.1 mL-1,200 mL
- User configurable maximum flow rate
- User configurable occlusion pressures
- Alarm before end
- Air bubble size detection
- KVO feature

To take advantage of this limited-time offer, contact your MEDRAD Sales Representative.

DF-112328-01